

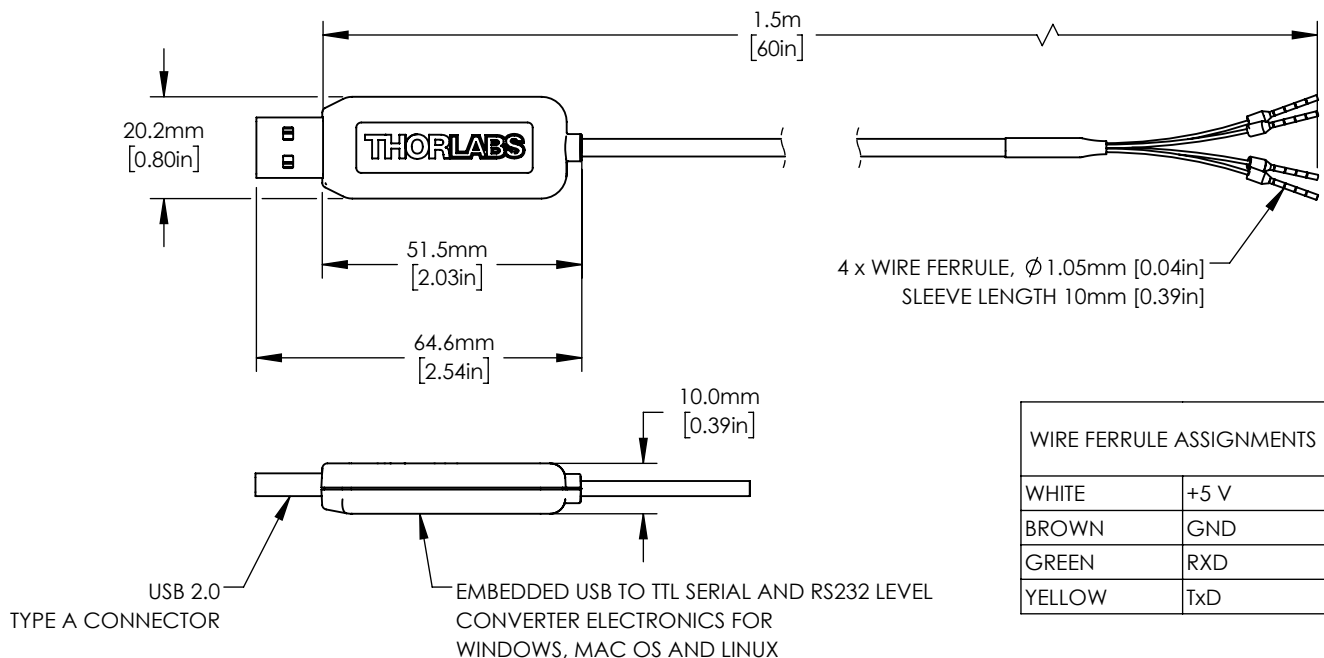


PM103E
PLUG WIRE FERRULES TO DFC 0.5/7 PUSH SPRING
CONNECTOR PROVIDED WITH PM103E

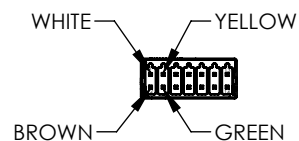


CONNECTOR PIN DIAGRAM	
1	+5 V (WHITE)
2	GND (BROWN)
3	TxD (YELLOW)
4	RxD (GREEN)
5-14	N.C.

PM101, PM102, PM103
PLUG WIRE FERRULES DIRECTLY TO THE D-SUB PINS

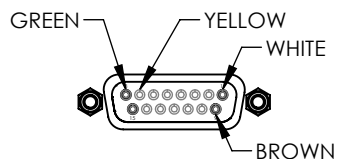

CONNECTOR PIN DIAGRAM	
1	+5 V (WHITE)
7	TxD (YELLOW)
8	RxD (GREEN)
9	GND (BROWN)
2-6, 10-15	N.C.

DEVICES WITH 9 PIN FEMALE RS232 CONNECTOR
LIKE PM101R AND PM5020
PLUG WIRE FERRULES DIRECTLY TO THE D-SUB PINS
ISOLATE THE NOT USED WHITE WIRE
PM101R MUST BE POWERED WITH USB CONNECTOR

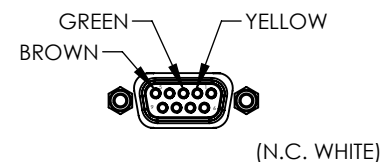

CONNECTOR PIN DIAGRAM	
2	TxD (YELLOW)
3	RxD (GREEN)
5	GND (BROWN)
1, 4, 6-9	N.C.

FOR INFORMATION ONLY
NOT FOR MANUFACTURING PURPOSES

DRAWING PROJECTION			 www.thorlabs.com USB A to Serial Flying Leads Adapter Cable	
NAME	DATE			
DRAWN	BHL	06/JUN/23	MATERIAL	
APPROVAL	MGO	12/JUN/23		
COPYRIGHT ©2023 BY THORLABS			N/A	
VALUES IN PARENTHESIS ARE CALCULATED AND MAY CONTAIN ROUND OFF ERRORS			ITEM #	APPROX WEIGHT
			CABRU15	50 g
			REV	A