DTSA02 - Spacer Plate for DTS25 and DTS50 Stages

For the continuing safety of the operators of the equipment and for the protection of the equipment itself, the operator should take note of the **Warnings, Cautions and Notes** throughout this leaflet, and, where applicable, in any associated literature.

1.1 Introduction

The DTSA02 spacer plate is used raise the deck height of the DTS25 and DTS50 stages to 30 mm.

..



Fig. 1.1 DTSA02 Spacer Plate, Mounted to a DTS50 Stage

1.2 Fitting .

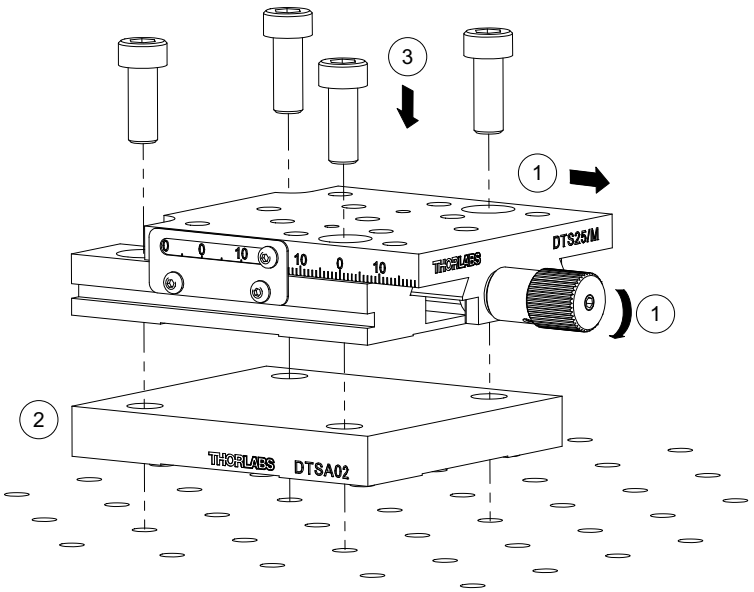


Fig. 1.2 Fitting the DTSA02 Spacer Plate

To fit the spacer plate, refer to Fig. 1.2 and proceed as follows:

1) Adjust the actuator to position the moving carriage such that the mounting holes in the base are clearly visible through the holes in the top plate.
2) Position the stage on the base plate, ensuring that the through holes in the spacer plate are clearly visable through the holes in the top plate.
3) Fit four bolts (M6 x 15 or 1/4-20 x 1/2”, not supplied) through the holes in the top platform and tighten to secure the stage to the work surface.

Thorlabs Worldwide Contacts

For technical support or sales inquiries, please visit us at www.thorlabs.com/contact for our most up-to-date contact information.



USA, Canada, and South America
Thorlabs, Inc.
sales@thorlabs.com
techsupport@thorlabs.com

Europe
Thorlabs GmbH
europe@thorlabs.com

France
Thorlabs SAS
sales.fr@thorlabs.com


Japan
Thorlabs Japan Inc.
sales@thorlabs.jp

UK and Ireland
Thorlabs Ltd.
sales@uk.thorlabs.com
techsupport.uk@thorlabs.com

Scandinavia
Thorlabs Sweden AB
scandinavia@thorlabs.com

Brazil
Thorlabs Vendas de Fotônicos Ltda.
brasil@thorlabs.com

China
Thorlabs China
chinasales@thorlabs.com

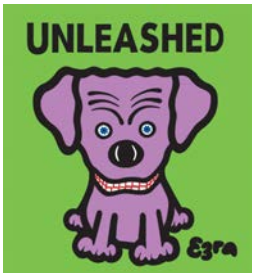

THORLABS
www.thorlabs.com