INF 0023 – Fitting Casters to a Breadboard Support Frame

This information sheet details how to fit the optional casters to a PFL series Breadboard Support Frame.

For the continuing safety of the operators of the equipment and for the protection of the equipment itself, the operator should take note of the **Warnings, Cautions and Notes** throughout this leaflet, and, where applicable, in any associated literature.

1 Parts List

Bolt qty (qty 4)

Washer (qty 4)

Nut (qty 4).

End Cap (qty 4)

Caster Assembly (qty 4)

2 Fitting the Caster Assemblies

Referring to Fig 2.1:

1) Fit a washer to the attachment bolt and pass the bolt through the holes in the lower support strut.
2) Using an open ended spanner, turn each mounting foot and adjust the height of the frame sufficiently to allow the caster assembly to be fitted.
3) Fit the caster assembly onto the bolt.
4) Fit the retaining nut and tighten.
5) Fit the plastic end cap to the head of the bolt.
6) Repeat items 1) to 6) for the remaining 3 casters.

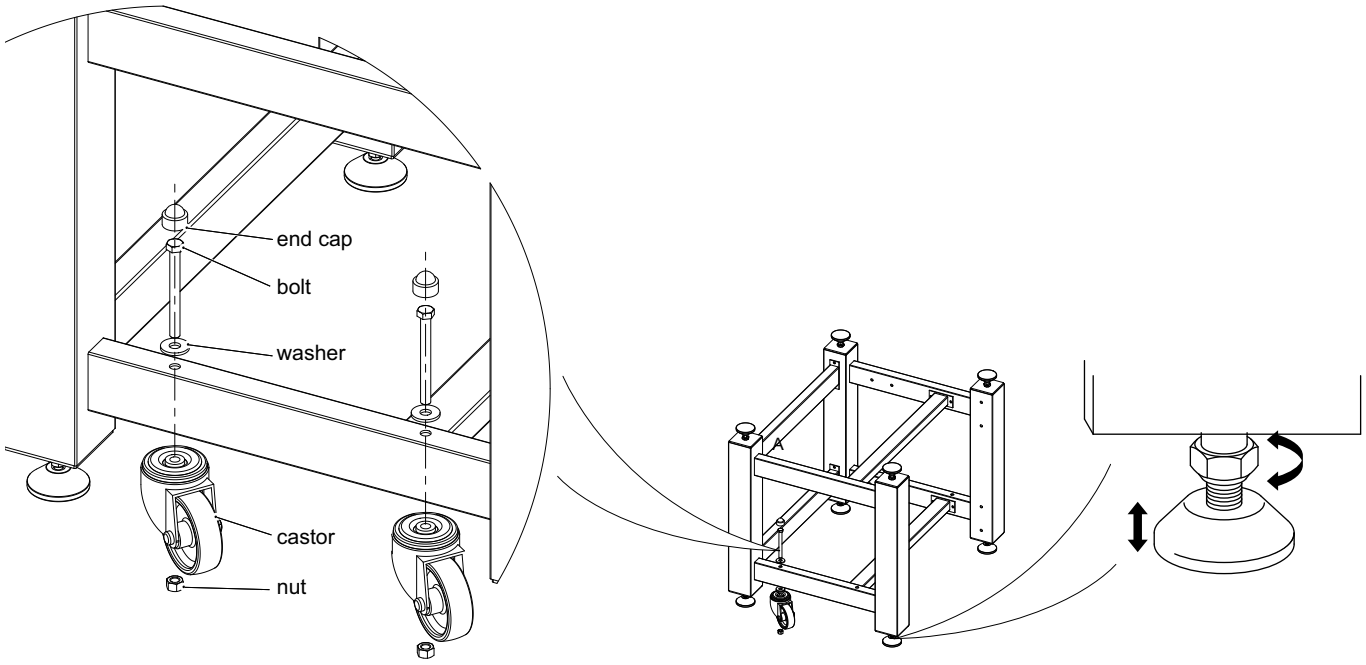


Fig. 2.1 Fitting a caster assembly

