Optics Cart Height Adjustment Pads - Fitting Instructions

This information sheet details instructions for fitting the height adjustment pads to a POC001 Optics Cart. For the continuing safety of the operators of the equipment and for the protection of the equipment itself, the operator should take note of the **Warnings, Cautions and Notes** throughout this leaflet, and, where applicable, in any associated literature.

1.1 Fitting Procedure

1) Unpack the cart and remove the packaging and protective wrapping.
2) Locate the four adjustment pads which are packed separately.
3) Position a pad on one of the height adjusters and locate the ball joint into the locating recess on the underside of the pad.
4) Press down to click the ball joint into the recess. Body weight should be sufficient but if the joint is tight, position a flat piece of wood on top of the pad and give a sharp tap with a mallet.
5) Repeat items (3) and (4) for the remaining height adjusters.

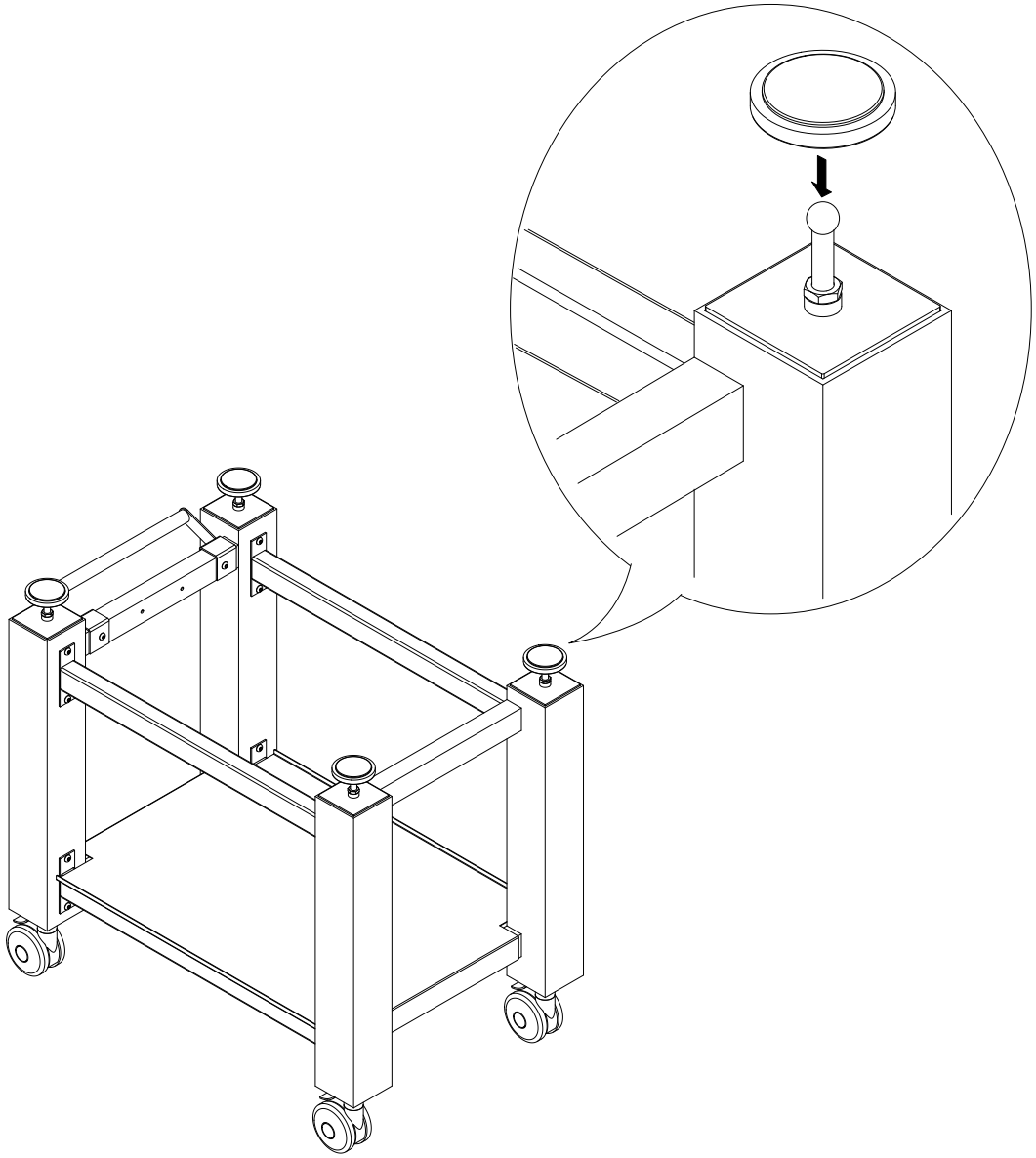


Fig. 1.1 Fitting the Height Adjustment Pads

Thorlabs Worldwide Contacts

For technical support or sales inquiries, please visit us at www.thorlabs.com/contact for our most up-to-date contact information.



USA, Canada, and South America
Thorlabs, Inc.
sales@thorlabs.com
techsupport@thorlabs.com

Europe
Thorlabs GmbH
europe@thorlabs.com

France
Thorlabs SAS
sales.fr@thorlabs.com


Japan
Thorlabs Japan Inc.
sales@thorlabs.jp

UK and Ireland
Thorlabs Ltd.
sales@uk.thorlabs.com
techsupport.uk@thorlabs.com

Scandinavia
Thorlabs Sweden AB
scandinavia@thorlabs.com

Brazil
Thorlabs Vendas de Fotônicos Ltda.
brasil@thorlabs.com

China
Thorlabs China
chinasales@thorlabs.com

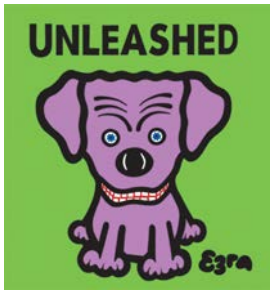

THORLABS

www.thorlabs.com