INF0136T – PWA100 and PWA110 Under Shelf - Fitting Instructions

For the continuing safety of the operators of the equipment and for the protection of the equipment itself, the operator should take note of the **Warnings, Cautions and Notes** throughout this leaflet, and, where applicable, in any associated literature.

1.1 Fitting the Mounting Brackets

Referring to Fig. 1.1:

1) Position the brackets on the underside of the shelf as shown below, such that the holes in the brackets align with the holes in the shelf, then using the eight countersunk screws and a suitable cross-head screwdriver, attach the brackets to the shelf. Do not over tighten or apply excessive force as this may cause the nut inserts to become detached..

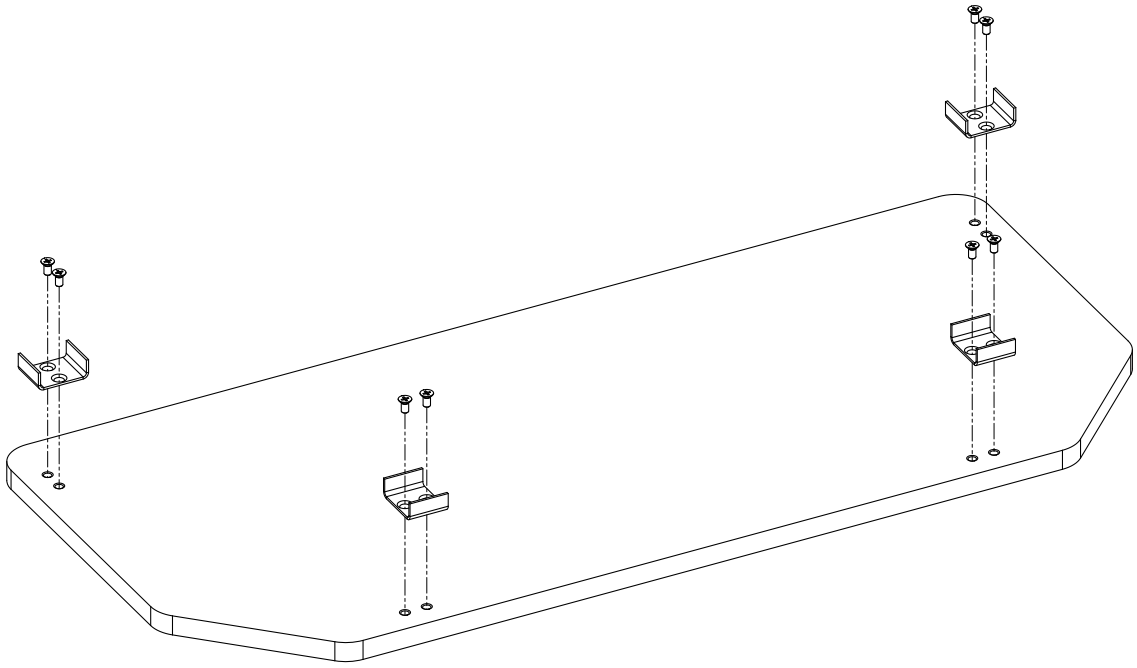


Fig. 1.1 Fitting the mounting brackets

1.2 Fitting the Shelf

Warning

When performing sub-item (1) below, good lifting practices should be observed and proper posture maintained. Do not lift loads in excess of 25kg per operator. If in doubt employ 2 persons.

1) Position the shelf on the lower rails frame, ensuring that the brackets sit correctly over the rails as shown below.

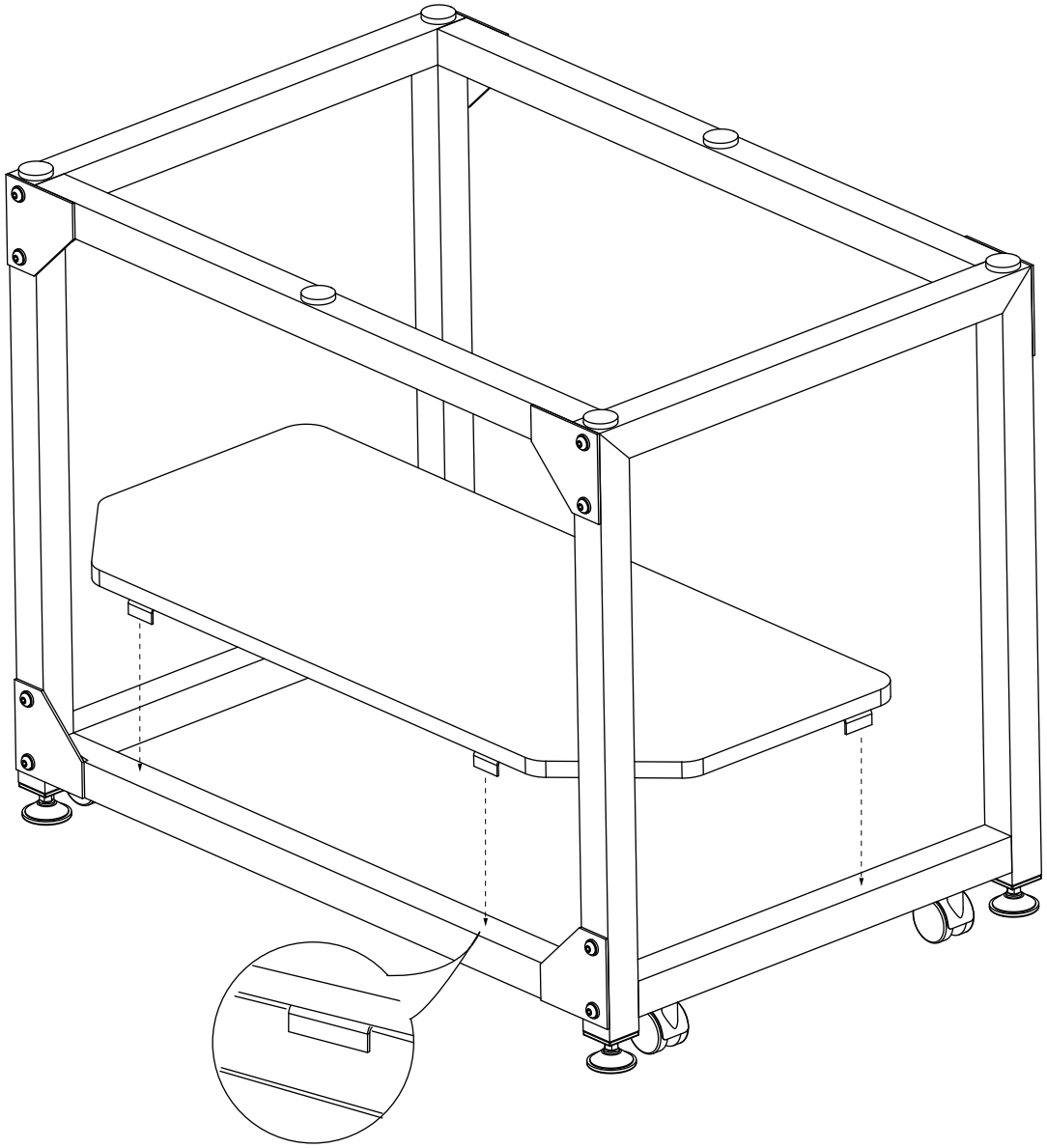


Fig. 1.2 Fitting the shelf

Thorlabs Worldwide Contacts

For technical support or sales inquiries, please visit us at www.thorlabs.com/contact for our most up-to-date contact information.



USA, Canada, and South America
Thorlabs, Inc.
sales@thorlabs.com
techsupport@thorlabs.com

Europe
Thorlabs GmbH
europe@thorlabs.com

France
Thorlabs SAS
sales.fr@thorlabs.com


Japan
Thorlabs Japan Inc.
sales@thorlabs.jp

UK and Ireland
Thorlabs Ltd.
sales@uk.thorlabs.com
techsupport.uk@thorlabs.com

Scandinavia
Thorlabs Sweden AB
scandinavia@thorlabs.com

Brazil
Thorlabs Vendas de Fotônicos Ltda.
brasil@thorlabs.com

China
Thorlabs China
chinasales@thorlabs.com

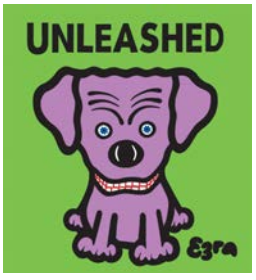

THORLABS
www.thorlabs.com